

VI ByPass Syringe™ Clear

- Designed to accommodate both long-neck and short-neck bottles
- A universal collection device
- Combines traditional products into a simple, intuitive closed system
- Designed to streamline blood collections, reduce false positive blood cultures and make blood collection easier from existing lines.

VI ByPass Syringe

19 gauge side channel permits 0-4 mL aspiration of waste volume

VI Velocity Reduction™

17 gauge cannula helps decrease device-related hemolysis

VI Sliding Fluid Seal™


Designed to improve sample quality by separating 99.9% of potential contaminants and helping reduce micro bubbles

Evaluated and Endorsed

by **ANTT**
Aseptic Non Touch Technique

(Aseptic Non Touch Technique)

Reduce contamination and ensure best practice

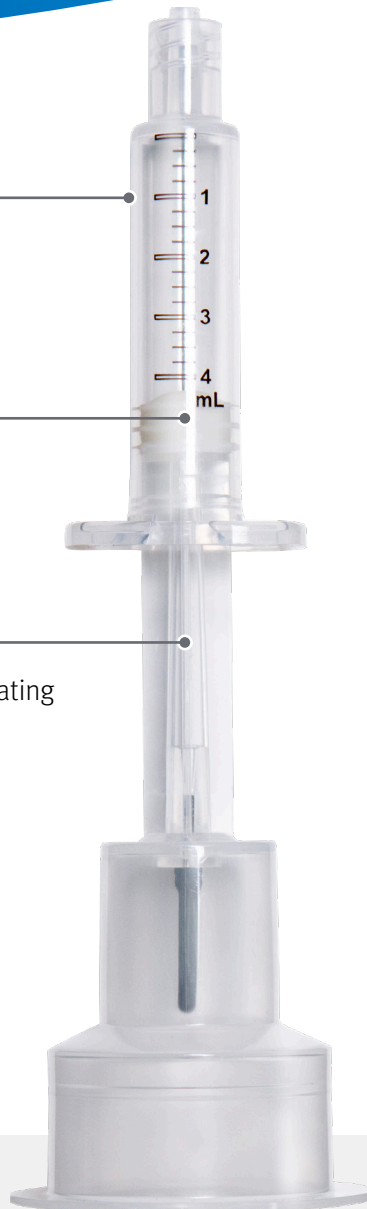

SKU# VI-101-WM-A

- Description: ByPass Syringe – Widemouth w/o Insert
- Unit/Box: 25
- Unit/Case: 100

Additional Details

- Direct Vacuum
- Diversion
 - Venipuncture
 - PIV insertion
- Vacuum source
 - Tubes
 - Bottles

Learn more by visiting www.vascularintegrity.com.